



## SCALP VEIN INFUSION SETS

# YOUR CHOICE DEFINES YOUR SAFETY NEED

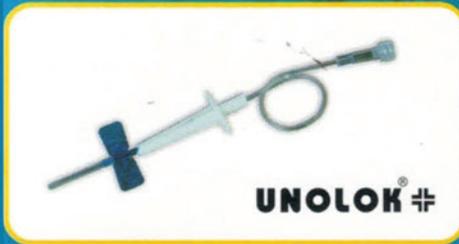

## Product Range

Needle Gauge	Colour Code	Needle Length (mm) & (inch)	Needle Dia. (mm)
18G	Pink	19 3/4"	1.20
19G	Cream	19 3/4"	1.10
20G	Yellow	19 3/4"	0.90
21G	Green	19 3/4"	0.80
22G	Black	19 3/4"	0.70
23G	Light Blue	19 3/4"	0.60
24G	Medium purple	19 3/4"	0.55
25G	Orange	19 3/4"	0.50
26G	Brown	19 3/4"	0.45

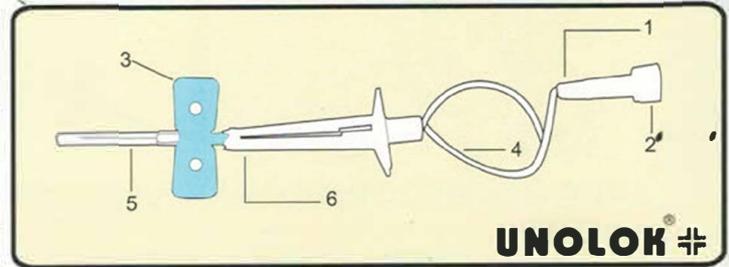

**5 NEEDLE** - The needle is thin walled which allows a larger internal diameter and is ideal for rapid infusion and dosage with high viscosity medications. The needle has a short bevel three facet point to minimize the chances of cross perforation of vein during insertion. Needles are siliconised to reduce friction & pain (for smooth penetration and withdrawal).

**6. PROTECTIVE SHIELD** in UNOLOK<sup>®</sup> † assists in ensuring there is no secondary needle sticks and no reuse of device

## Components

**1. Luer Adaptor** - Luer adaptors are transparent which allow visual monitoring of the flow of the medicaments. Luer adaptors are having a 6% luer taper as per ISO-594 which provides a leak free mating with an infusion set or a syringe of any make and size.

**2. LUER LOCK CONNECTOR:** Presence of rigid PVC luer lock connector in UNOLOK Scalpvein sets provides additional secure fitment and prevents any accidental leakage.

**3. WINGS** - The soft wings conform well to the skin contours and make immobilization simple.

**4. LINK TUBE** - The link tube is translucent for easy observation of the contents.

## PACKING

Individually packed in transparent thermoformed peelable blister pack.

	UNIT PACK	SHELF CARTON	SHIPPING CORRUGATED BOX
UNOLOK <sup>+</sup>	Ten blister of UNOLOK <sup>+</sup> S.V. sets from a strip each unit is easily separable from the strip along with perforated slitting	100 units of UNOLOK <sup>+</sup> S.V. Sets (16 strips of six units and one strip of four units each are packed in a printed duplex board box.)	10 shelf cartons are packed in a 5 ply corrugated cartons. Pack Size: 360x330x495mm Gross Wt.: 6.0 Kg. (approx)