




# Cathula<sup>®</sup>†

Safety I.V. Cannula with *tipclip*<sup>™</sup>  
*Protecting the protectors*

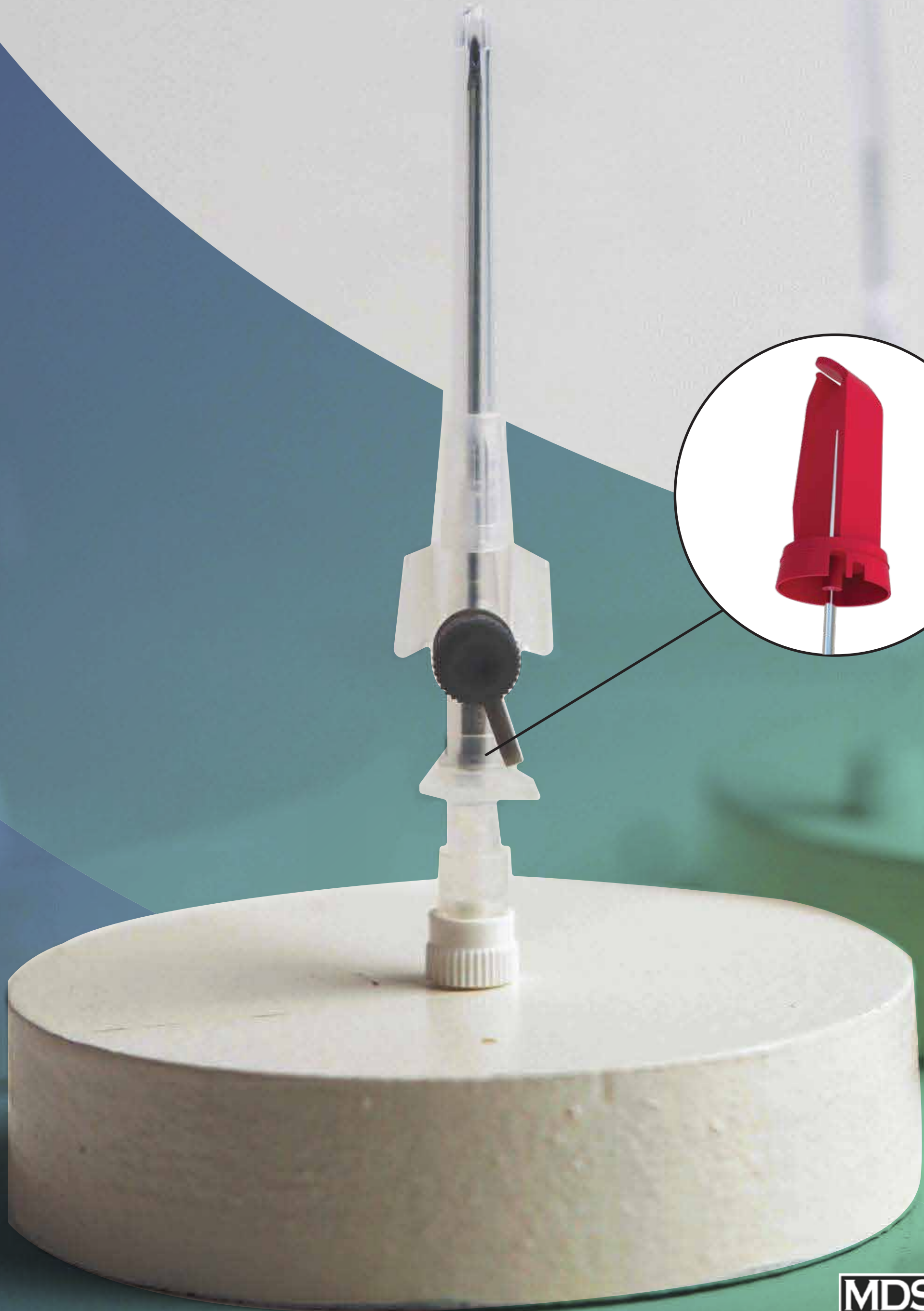

HINDUSTAN SYRINGES  
& MEDICAL DEVICES LTD.  
ADDRESS ISO 13485 : 2016

*As part of your job and your passion, you care for hundreds of patients every day. Doctor, nurses and medical staff are depended upon in time of need. As a health professional you can now depend upon someone as well.*

## HEALTHCARE PROFESSIONALS ARE AT HIGH RISK

As per WHO, among the 35 million health professionals worldwide, about 3 million receive percutaneous exposures to blood borne pathogens each year, 90% of these infections occur in developing countries.

## WHENEVER YOU ARE ON DUTY, CATHULA + IS ALSO ON DUTY

A protection device **tipclip™** encapsulates the tip of the used needle when introducer needle is extracted from the Catheter; Post Cannulation. **tipclip™** prevents accidental needle stick injury provided that general effective and safe working procedures and precautions are maintained during its use and disposal

## AUTO PROTECTION

The hassel-free and auto protection technology of Cathula+ I.V.Cannula is designed to provide 'Involuntary Activated Safety'-an intuitive mechanism that requires no additional training of the health staff and strives to make I.V. insertion and introducer needle extraction completely safe.

### INTRODUCER NEEDLE

Siliconised electrofinished, back cut bevelled stainless steel needle.

Instant Flash Back Notch on 24G and 26G Cannula.

**tipclip™**

Sharp Injury  
Prevention Clip.



### NEEDLE HUB

Allows blood flash back visualisation.

Thumb rest for proper grip during insertion.

### LUER LOCK CAP

Closure to wing housing permits patient mobility when infusion temporarily discontinued.

### BLOOD STOPPER

Designed to expedite blood flash back for enabling correct and efficient cannulation.

### BLISTER PACKAGE

Cannula is packed in rigid blister and sealed in peel open type medical grade paper.

### CATHETER

Available in FEP radio detectable with 4 strip

### INJECTION PORT CAP

It facilitates extra medication while preventing its back flow.



**MDSAP**  
MEDICAL DEVICE SINGLE AUDIT PROGRAM



Scan to  
Visit Product Video



Scan to  
Visit Product Page

COMPONENTS