



**DISPO VAN**<sup>®</sup>

**SINGLE USE SYRINGES**

---

Most Trusted Brand of India






## STERILE SINGLE USE SYRINGES

Dispovan single use hypodermic syringes is the largest selling single-use syringes in the Indian market. This sterile hypodermic syringe for single use is manufactured as per IS:10258/ISO:7886-1 (Sterile Hypodermic Syringe for Single use) for normal use. Furthermore, this 3-piece syringe has a special inert Piston which minimizes chances of leakage. It also has a shelf life of 5 years.

## COMPONENTS

### **BARRELS AND PLUNGERS:**

Made of non-toxic, medical-grade polypropylene compatible with the medication for short term contact drug delivery application. It comes with a luer mount nozzle with a 6% luer taper as per ISO - 80369-7

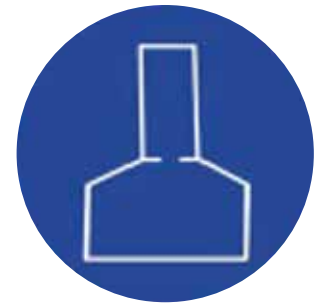

**PACKING:** DISPO VAN single use syringes are individually Ribbon packed in a robust double-laminated plastic film. Ribbon packaging is superior to blister packaging due to higher puncture and water resistance, enabling extended shelf-life under adverse conditions. But for enabling optimum container loadability it is also available in transparent thermoformed peelable Blister pack with medical grade paper.

### **GASKETS:**

Made of natural rubber or thermoplastic elastomer which is chemically inert and compatible with medications for short term contact, drug delivery applications and for ensuring smooth movement of plunger.

# THE STERILE HYPODERMIC SYRINGES FOR SINGLE USE ARE PREFERRED FOR THE FOLLOWING REASONS:

- Manufactured on highly automated machines with vision inspection.
- Prevent the spread of blood-borne infections and illnesses.
- Mass-produced, allowing greater accessibility and easy availability, meeting any urgent requirements.
- Do not require cleaning or sharpening, leading to low-maintenance costs.
- Affordable and widely available.
- Cannula made with premium stainless steel and engineered for a conical and sharp bevel ensuring comfortable insertion.

CAPACITY	TYPE OF NOZZLE	STANDARD NEEDLE DIMENSIONS
1 ML. * (T.B)	Luer mount/ luer lock with central nozzle	26G X 1/2 <sup>Ø</sup> *0.45 x 13
2 ML. **	Luer mount/ luer lock with central nozzle	24G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 22G X 1 <sup>Ø</sup> *0.55 x 25   *0.60 x 25   *0.60 x 38   *0.70 x 25 22G x 1 1/4 <sup>Ø</sup> 22G x 1 1/2 <sup>Ø</sup> 21G x 1 <sup>Ø</sup> 21G x 1 1/2 <sup>Ø</sup> *0.70 x 32   *0.70 x 38   *0.80 x 25   *0.80 x 38
3ML	Luer mount/ luer lock with central nozzle	24G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 22G X 1 <sup>Ø</sup> *0.55 x 25   *0.60 x 25   *0.60 x 38   *0.70 x 25 22G x 1 1/4 <sup>Ø</sup> 22G x 1 1/2 <sup>Ø</sup> 21G x 1 <sup>Ø</sup> 21G x 1 1/2 <sup>Ø</sup> *0.70 x 32   *0.70 x 38   *0.80 x 25   *0.80 x 38
5ML	Luer mount/ luer lock with central nozzle	24G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 23G X 1 <sup>Ø</sup> 22G X 1 <sup>Ø</sup> *0.55 x 25   *0.60 x 25   *0.60 x 38   *0.70 x 25 22G x 1 1/4 <sup>Ø</sup> 22G x 1 1/2 <sup>Ø</sup> 21G x 1 <sup>Ø</sup> 21G x 1 1/2 <sup>Ø</sup> *0.70 x 32   *0.70 x 38   *0.80 x 25   *0.80 x 38
10ML**	Luer mount with eccentric nozzle    Luer Lock with central nozzle	22G X 1 1/4 <sup>Ø</sup> 22G X 1 1/2 <sup>Ø</sup> 0.70 x 32    0.70 x 38 *21G x 1 1/2 <sup>Ø</sup> 20G x 1 1/2 <sup>Ø</sup> 0.80 x 38    0.90 x 38
20ML**	Luer mount with eccentric nozzle    Luer Lock with central nozzle	22G X 1 1/4 <sup>Ø</sup> 22G X 1 1/2 <sup>Ø</sup> 0.70 x 32    0.70 x 38 *21G x 1 1/2 <sup>Ø</sup> 20G x 1 1/2 <sup>Ø</sup> 0.80 x 38    0.90 x 38
50ML	Luer mount with eccentric nozzle    Luer Lock with central nozzle	Without Needle



Conform to : ISO:7886-1/IS: I 0258 standards (Sterile hypodermic syringe for single use).  
Before use check the integrity of the sealing of the pack.  
Do not use if previously opened or damage.



Scan to  
Visit Product Page