

SCALP VAN




SCALP VEIN INFUSION SETS

**YOUR CHOICE
DEFINES YOUR
COMFORT NEEDS**



SCALP VAN

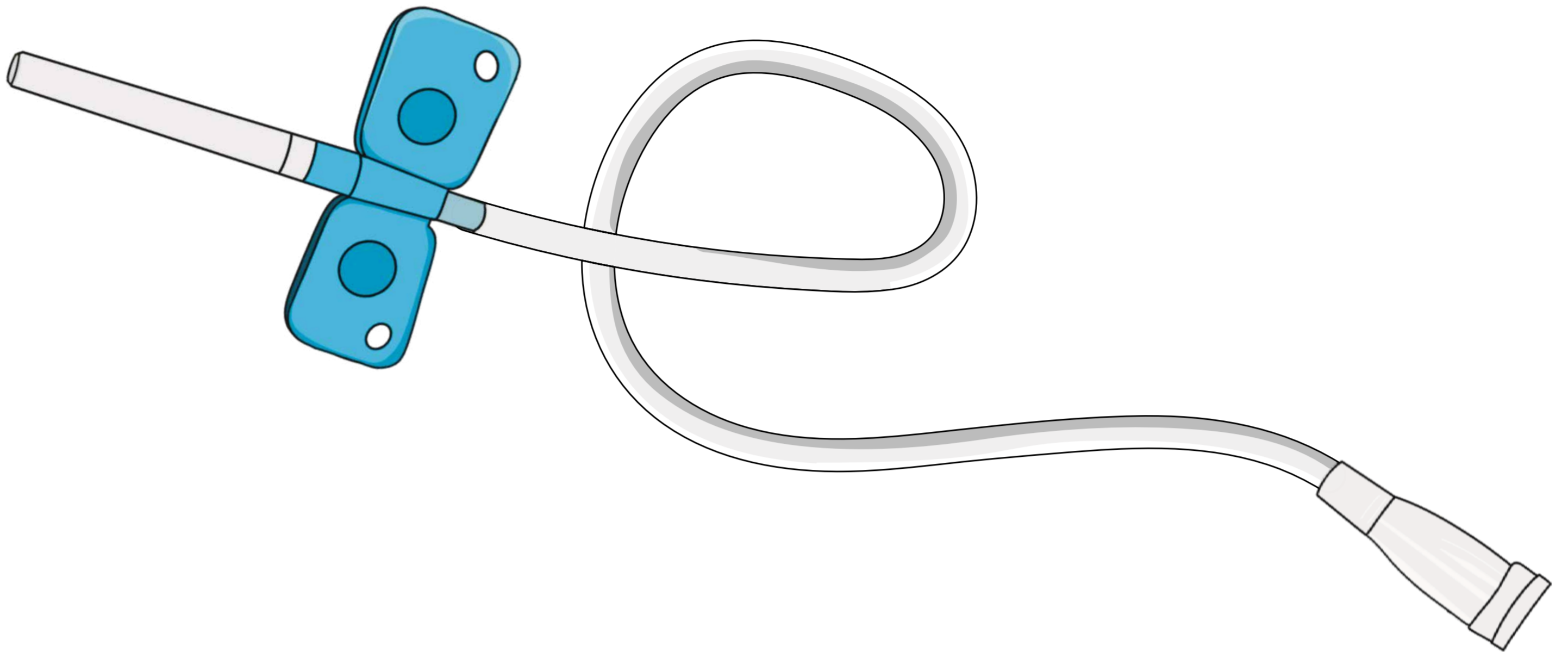

COMPONENTS

LUER ADAPTOR

- Luer adaptors are transparent which allow visual monitoring of the flow of the medicaments.
- Luer adaptors have a 6% luer taper as per ISO-594 which provides a leak free mating with an infusion set or a syringe.

WINGS

- The soft wings conform well to the skin contours and make immobilisation simple.

KINK TUBE

- The kink free link tube is translucent for easy observation of fluids for instant flashback to enable quick cannulation with minimal risk of cross perforation of vein.

NEEDLE

- The needle is thin walled which allows a larger internal diameter and is ideal for rapid infusion.
- The needle has a short bevel to minimise the chances of cross perforation of vein during insertion.
- Needles are siliconised to reduce friction & pain (for smooth penetration and withdrawal)

PRODUCT RANGE

Needle Guage	Colour Code	Needle Length		Needle Diameter (mm)
		(mm)	(Inches)	
18G	Pink	19	3/4 "	1.20
19G	Cream	19	3/4"	1.10
20G	Yellow	19	3/4"	0.90
21G	Green	19	3/4"	0.80
22G	Black	19	3/4"	0.70
23G	Light Blue	19	3/4"	0.60
24G	Red	19	3/4"	0.55
25G	Orange	19	3/4"	0.50
26G	Brown	19	3/4"	0.45

PACKING

Individually packed in transparent thermoformed peelable blister pack

UNIT PACK	SHELF CARTON	SHIPPING CORRUGATED BOX
Two blister of SCALP VAN S.V sets form a strip, each unit is easily separable from the strip along with perforated slitting	50 units of SCALP VAN S.V sets (25 strips of two units each are packed in a printed duplex board box)	40 shelf cartons are packed in a 5 ply corrugated cartons. Pack Size 570x465x320 Gross Wt: 11.5 Kg (approx)

WARNING

Do not use if the pack is previously opened or damaged
Do not attempt to re-steralise and use
Infusion sets should not be recycled and should be disposed off preferably in sharps container.

PRECAUTIONS

Before use check the integrity of the sealing of the pack. Tape the wings securely after insertion.

