


# UNOLOK<sup>®</sup>

## SCALP VEIN INFUSION SETS

YOUR CHOICE FOR SECURE FITMENT

### LUER LOCK CONNECTOR

A rigid PVC luer lock connector in UNOLOK Scalp vein winged infusion sets provides additional security of fitment to prevent any accidental leakage.

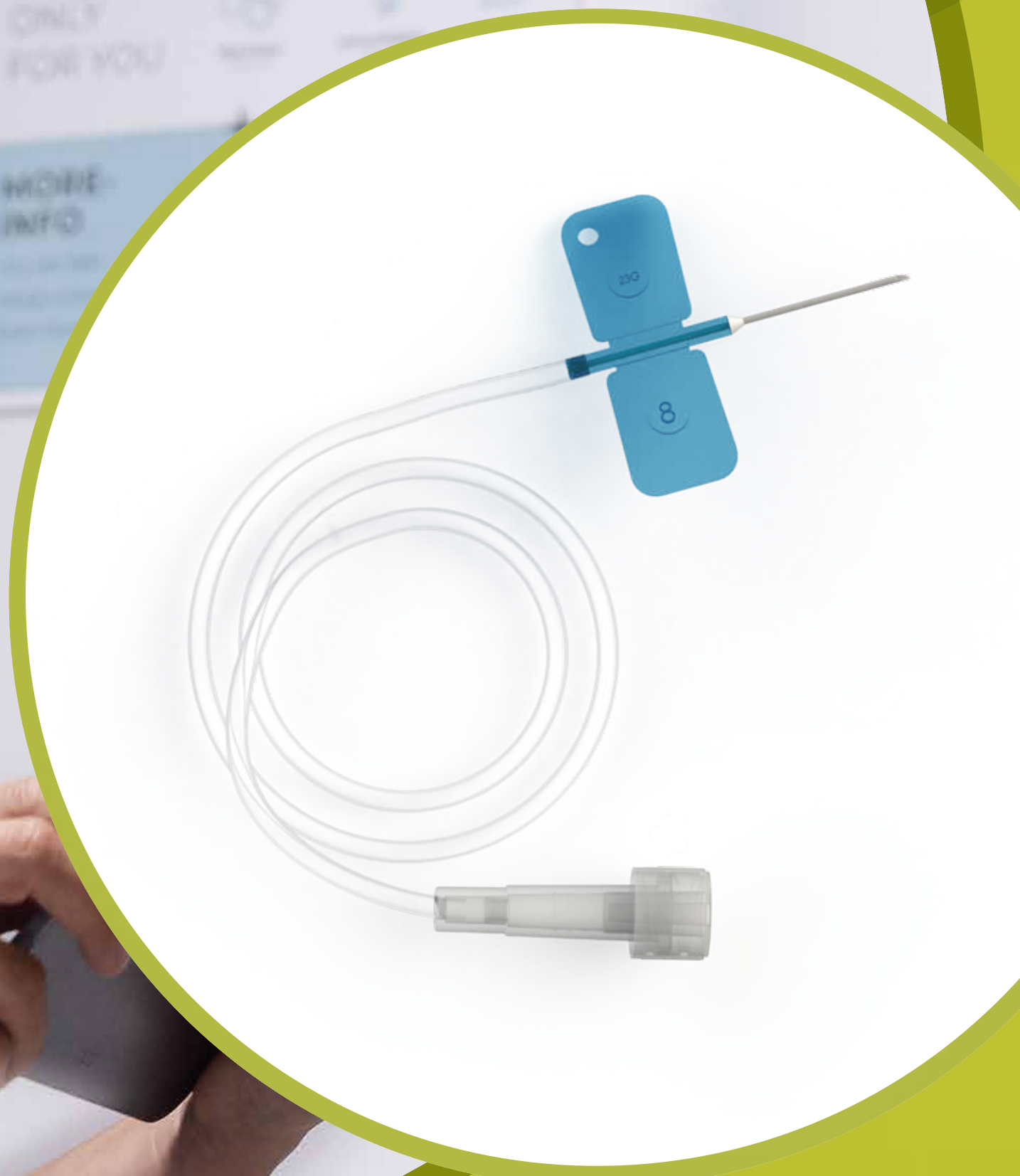



# PRODUCT RANGE

Needle Guage	Colour Code	Needle Length		Needle Diameter (mm)
		(mm)	(Inches)	
18G	Pink	19	3/4 "	1.20
19G	Cream	19	3/4"	1.10
20G	Yellow	19	3/4"	0.90
21G	Green	19	3/4"	0.80
22G	Black	19	3/4"	0.70
23G	Light Blue	19	3/4"	0.60
24G	Red	19	3/4"	0.55
25G	Orange	19	3/4"	0.50
26G	Brown	19	3/4"	0.45

# COMPONENTS

## LUER LOCK ADAPTOR

Luer adaptors are transparent which allows for visual monitoring of the flow of the medicaments. Luer Lock adaptors have a 6% luer taper as per ISO-594 which provides for a secure slip free mating with an infusion set or a syringe of any make and size to prevent accidental leakages.

## WINGS

The soft flexible wings conform well to the skin countours and make immobalisation simple.

## LINK TUBE

The kink free link tube is of high quality PVC material is translucent for easy observation of contents and enable instant flash back to aid cannulation.


## NEEDLE

The needle is thin walled which allows for a larger internal diameter and is ideal for rapid infusion and dosage with high viscosity medications. The needle has a short bevel and three facet point to minimise the chances of cross perforation of vein during insertion. Needles are siliconized to reduce friction & pain (for smooth penetration and withdrawal) to facilitate a comfortable insertion

# PACKING

Individually packed in transparent thermoformed peelable blister pack using superior quality medical grade paper

UNIT PACK	SHELF CARTON	SHIPPING CORRUGATED BOX
Ten blister of UNOLOK Scalp Vein sets from a strip each unit is easily separable from the strip along with perforated slitting.	100 units of UNOLOK Scalp Vein sets (10 strips of ten units each are packed in a printed duplex board box)	10 shelf cartons are packed in a 5 ply corrugated cartons.  Size: 355x3 10x395mm Gross Wt: 5.0 Kg (approx)

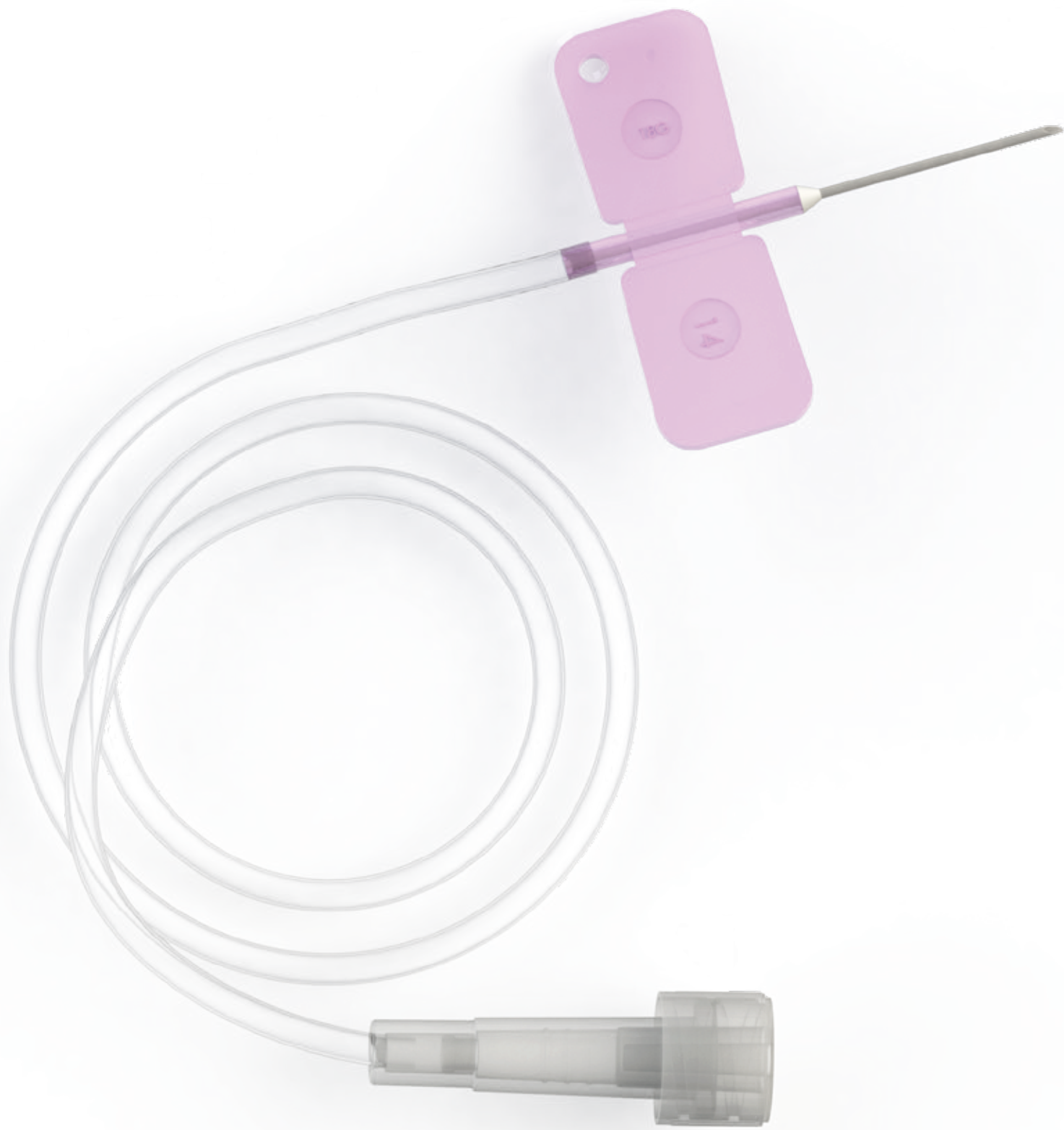

## WARNING

- Do not use if the pack is previously opened or damaged
- Do not attempt to re-sterilize and use
- Infusion sets should not be recycled and should be disposed off preferable in sharps containers

## PRECAUTIONS

Before use check the integrity of the sealing of the pack. Tape the wings securely after insertion.

