

Cathy[®] 

NEO with **INSTANT-FLSH BK**



Getting it Right 1st Time, Every Time

WE CARE FOR YOUR LITTLE ONES



As part of your job and your passion, you care for hundreds of patients every day. Doctor, nurses and medical staff are depended upon in times of need. As a health professional you can now depend on someone as well.

Whenever you are on duty, CATHY is also on duty

A protection device **SIPCLIP™** encapsulates the tip of the used needle when introducer needle is extracted from the catheter; Post Cannulation. **SIPCLIP™** prevents accidental needle stick injury provided that general effective and safe working procedures and precautions are maintained during its use and disposal.

Healthcare professionals are at high risk

As per WHO, among the 35 million health professionals worldwide, about 3 million receive percutaneous exposures to blood borne pathogens each year, 90% of these infections occur in developing countries.

Auto protection

The hassel-free and auto protection technology of CATHY I.V Cannula is designed to provide "Involuntary Activated Safety" – an intuitive mechanism that requires no additional training of the health staff and strives to make I.V insertion and introducer needle extraction completely safe.

CATHY NEO with "INSTANT FLSH BK" Eases Cannulation, Reduces Anxiety

For neonates and children, in 24G & 26G I.V Cannula, a notch has been provided in the needle for "INSTANT FLSH BK" (instant flashback visualisation) at the tip of the catheter. It assists in achieving successful Cannulation in 1st attempt of vein entry. It is specially helpful for neonates with difficult to access veins thus reducing trauma and anxiety – both for you and your patient.

NEEDLE

Siliconised electrofinished, back cut bevelled stainless steel needle
Notch on 24G and 26G Cannula

SIPCLIP™
Sharp Injury
Prevention Clip.




LUER LOCK CAP

Closure to wing housing permits patient mobility when infusion temporarily discontinued.

BLOOD STOPPER

The stopper has fine air vent to expedite blood flash back for enabling correct cannulation.

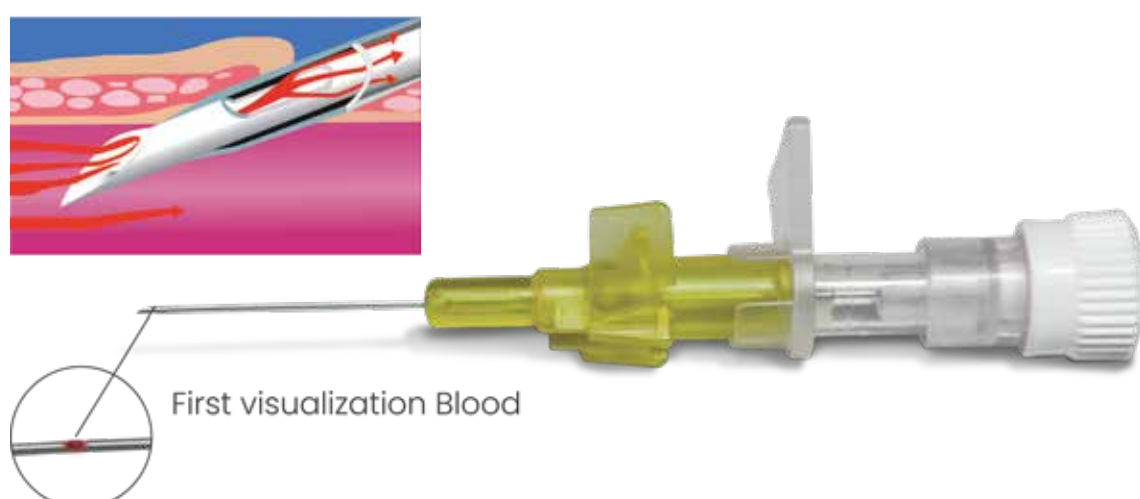

First visualization Blood

NEEDLE HUB

Allows blood flash back visualisation.
Thumb rest for proper grip during insertion.

CATHETER

Available in VOLEX™ polyuerthane radio detectable.

INJECTION PORT CAP

It facilitates extra medication and also prevents its back flow.

BLISTER PACKAGE

Cannula is packed in rigid blister and sealed in peel open type medical grade tyvek.

Colour Code	Guage	Catheter Ext. DiaxLength(mm)	Water flow-rate (ml/min)
Yellow	24G	0.70 x 19	16/12
Violet	26G	0.6 x 19	12/15



Scan to
Visit Product Video



Scan to
Visit Product Page