


Cathy[®] NEO

INSTANT-FLSH BK



Getting it Right 1st Time, Every Time

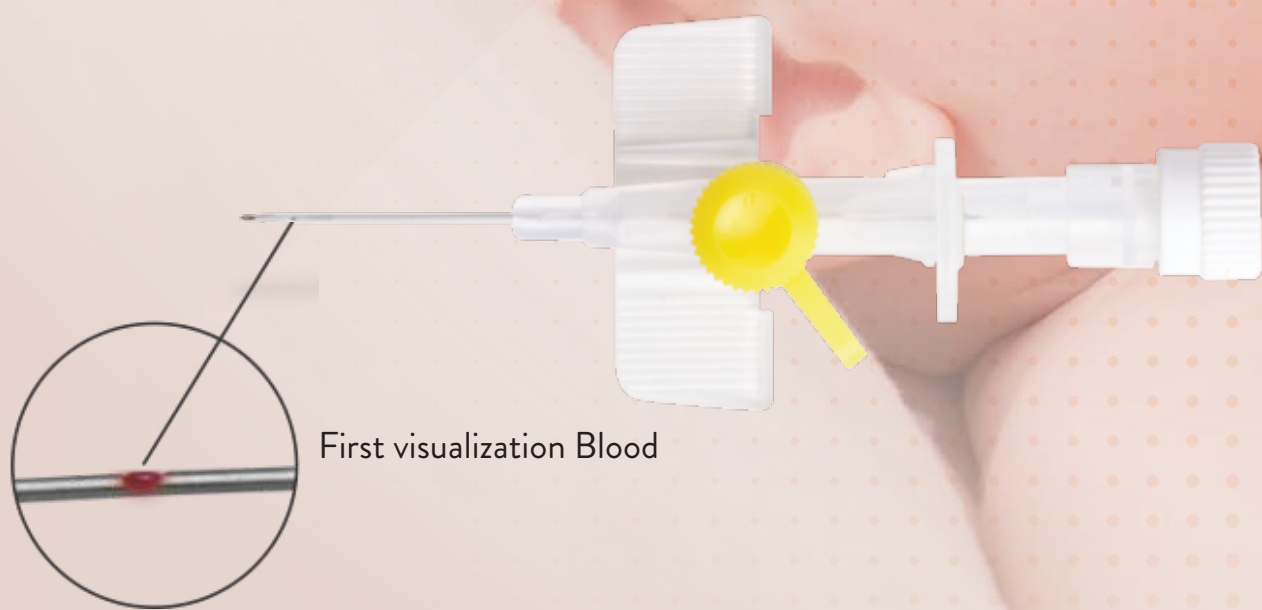

WE CARE FOR YOUR LITTLE ONES



CATHY with "INSTANT FLSH BK" Eases Cannulation, Reduces Anxiety

For neonates and children, in 24G & 26G I.V Cannula, a notch has been provided in the needle for "INSTANT FLSH BK" (instant flashback visualisation) at the tip of the catheter. It assists in achieving successful Cannulation in 1st attempt of vein entry. It is specially helpful for neonates with difficult to access veins thus reducing trauma and anxiety – both for you and your patient.

NEEDLE

Siliconised electrofinished, back cut bevelled stainless steel needle
Notch on 24G and 26G Cannula

LUER LOCK CAP

Closure to wing housing permits patient mobility when infusion temporarily discontinued.

BLOOD STOPPER

The stopper has fine air vent to expedite blood flash back for enabling correct cannulation.

NEEDLE HUB

Allows blood flash back visualisation.
Thumb rest for proper grip during insertion.

CATHETER



Available in VOLEX™ Polyurethane radio detectable 4 strips.

INJECTION PORT CAP

It facilitates extra medication and also prevents its back flow.

BLISTER PACKAGE

Cannula is packed in rigid blister and sealed in peel open type medical grade tyvek.

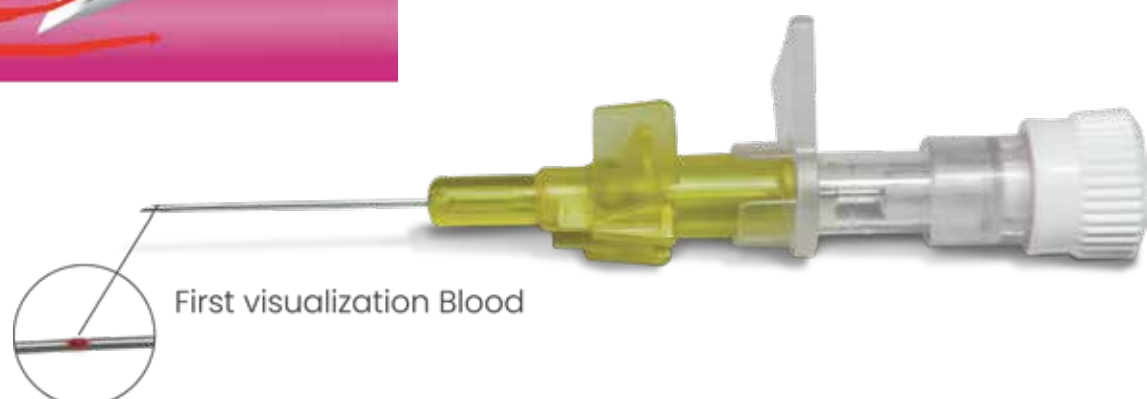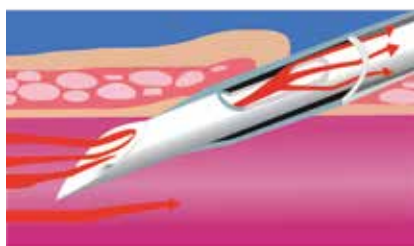

Colour Code	Gauge	Catheter Ext. Dia x Length(mm)	Water flow-rate (ml/min)
Yellow	24G	0.7 x 19	22
Violet	26G	0.6 x 19	15



Scan to
Visit Product Video



Scan to
Visit Product Page

Revision No:(Jun. 2024/Rev.-01)



HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.

174/25, Ballabgarh, Faridabad (Haryana), Pincode-121004. India | www.hmdhealthcare.com | Info@hmdhealthcare.com |