


SCALP VAN



SCALP VEIN INFUSION SETS

YOUR CHOICE
DEFINES YOUR
COMFORT NEEDS

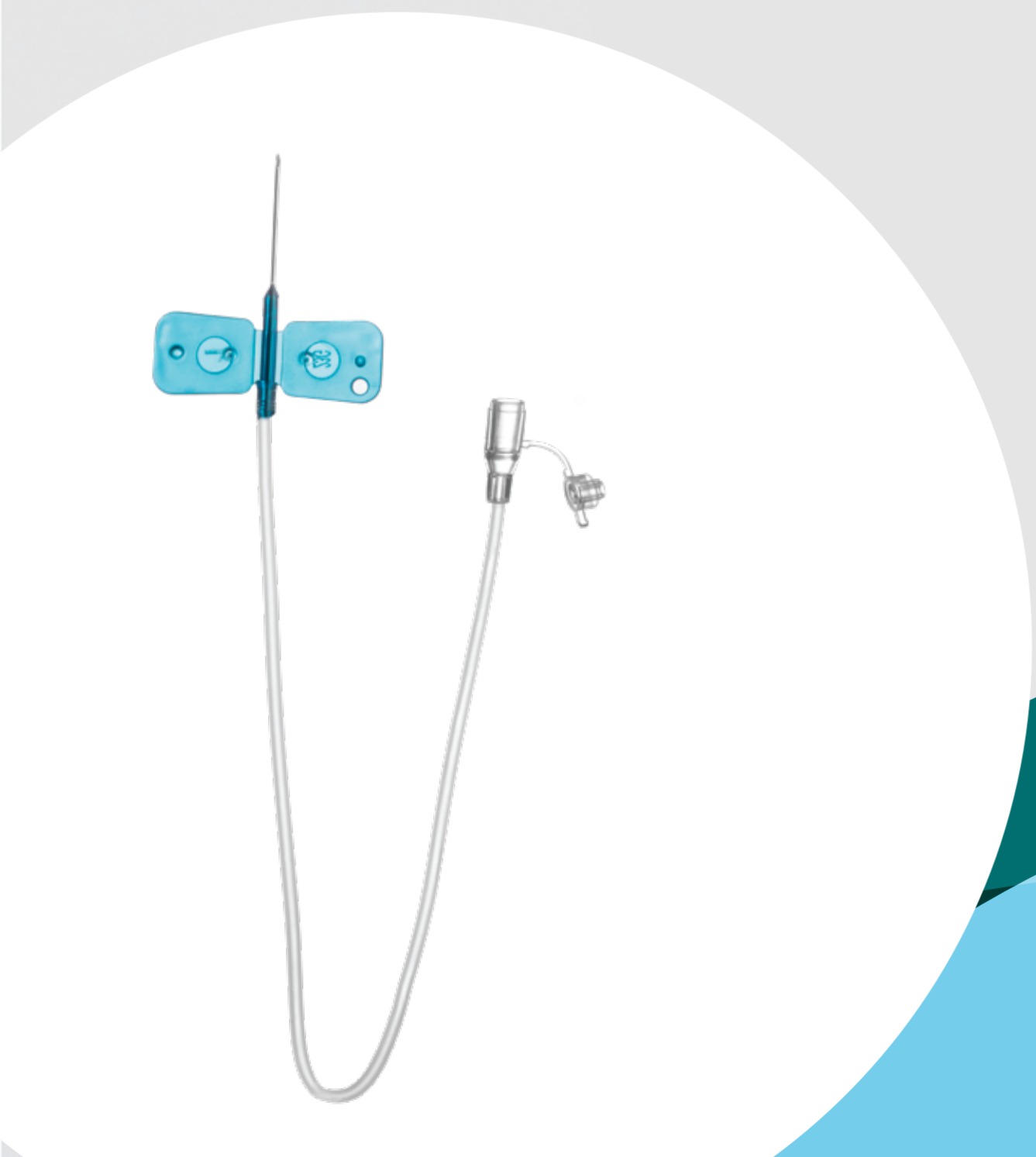

FEATURES

- Luer mount adapter with a 6% luer taper, are transparent which allows for visual monitoring of the flow of the medicaments & secure slip free mating prevent accidental leakages.
- The soft flexible wings conform well to the skin contours and make immobilization simple.
- The kink free link tube is of high quality PVC material is translucent for easy observation of contents and enable instant flash back to aid cannulation.
- The needle is thin walled which allows for a larger internal diameter and is ideal for rapid infusion and dosage with high viscosity medications.
- Short bevel and three facet needles are siliconized to reduce friction & pain to facilitate a comfortable insertion

PRODUCT RANGE

Needle Gauge	Colour Code	Needle Length		Needle Diameter (mm)
		(mm)	(Inches)	
18G	Pink	19	3/4 "	1.20
19G	Cream	19	3/4"	1.10
20G	Yellow	19	3/4"	0.90
21G	Green	19	3/4"	0.80
22G	Black	19	3/4"	0.70
23G	Light Blue	19	3/4"	0.60
24G	Purple	19	3/4"	0.55
25G	Orange	19	3/4"	0.50
26G	Brown	19	3/4"	0.45
27G	Grey	19	3/4"	0.40

PACKING

Individually packed in transparent thermoformed peelable blister pack

UNIT PACK	SHELF CARTON	SHIPPING CORRUGATED BOX
Two blister of SCALP VAN S.V. sets form a strip, each unit is easily separable from the strip along with perforated slitting	50 units of SCALP VAN S.V sets (25 strips of two units each are packed in a printed duplex board box)	40 shelf cartons are packed in a 5 ply corrugated cartons. Pack Size: 570x465x320 Gross Wt: 11.5 Kg (approx)

WARNING

Do not use if the pack is previously opened or damaged
Do not attempt to re-sterilise and use
Infusion sets should not be recycled and should be disposed off preferable in sharps containers

PRECAUTIONS

Before use check the integrity of the sealing of the pack. Tape the wings securely after insertion.



Scan to Visit Product
Revision No: Jun. 2024/Rev-01



HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.

174/25, Ballabgarh, Faridabad (Haryana), Pincode-121004. India | www.hmdhealthcare.com | Info@hmdhealthcare.com |