


Cathula[®] NEO



I.V. Cannula with “Instant Flush BK”

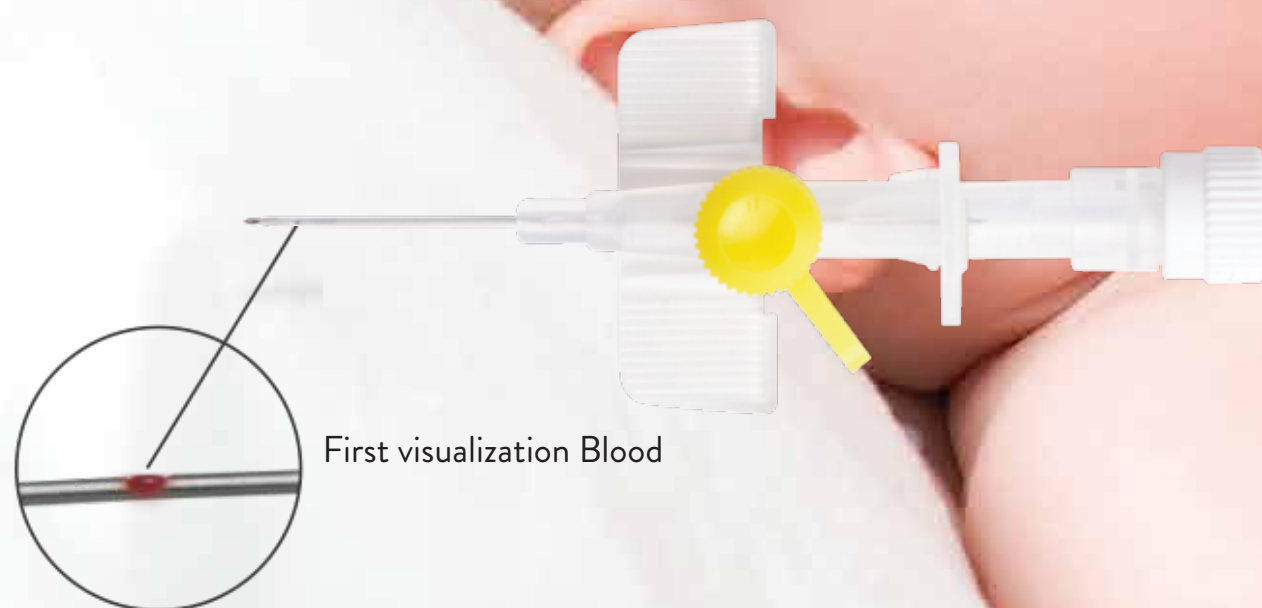

First visualization Blood

CATHULA NEO


To keep the little ones safe.



As part of your job and your passion, you care for hundreds of patients every day. Doctor, nurses and medical staff are depended upon in time of need. As a health professional you can now depend upon someone as well.

CATHULA· NEO WITH "INSTANT FLSH BK"

For neonates and children, in 24G & 26G IV. Cannula, a notch has been provided in the needle for "INSTANT FLSH BK" (instant flashback visualization) at the tip of the catheter to facilitate successful cannulation in 1st attempt of vein entry.

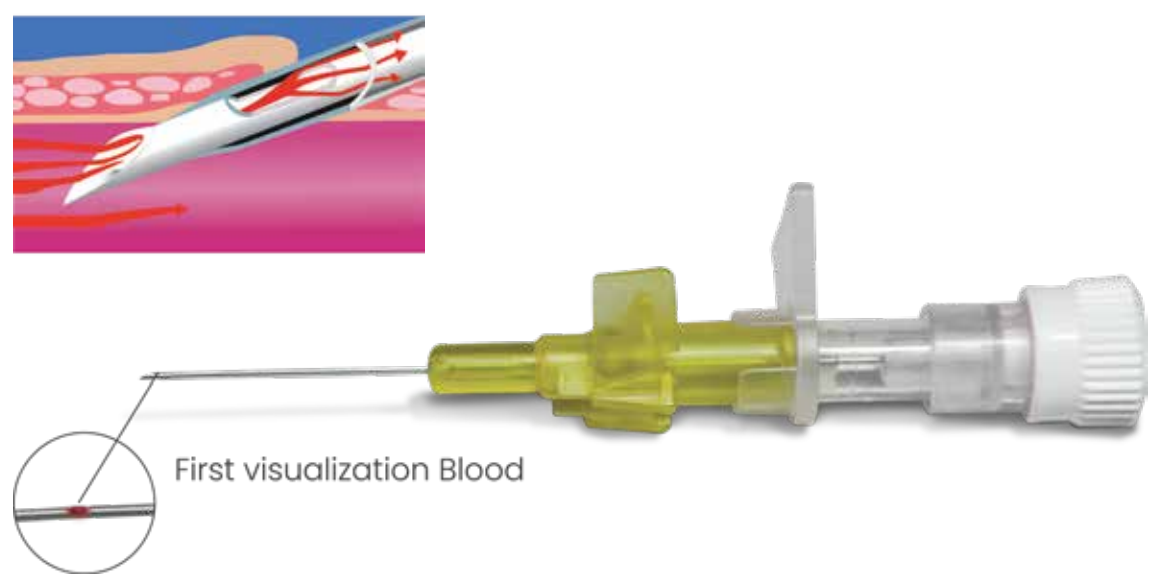

Colour Code	Gauge	Catheter Ext. Diameter x Length MM	Water flow-rate ml/min
Yellow	24G	0.7 x 19	22
Violet	26G	0.6 x 19	15

Reference: ISO 10555-1 & ISO: 10555-5

NEEDLE

Siliconised electrofinished, back cut bevelled stainless steel needle.
Instant Flash Back Notch on 24G and 26G Cannula.

LUER LOCK CAP

Closure to wing housing permits patient mobility when infusion temporarily discontinued.

BLOOD STOPPER

The stopper has fine air vent to expedite blood flash back for enabling correct cannulation.

NEEDLE HUB

Allows blood flash back visualisation.
Thumb rest for proper grip during insertion.

CATHETER

Available in FEP with radio detectable 4 strips

INJECTION PORT CAP

It facilitates extra medication by preventing its back flow.

BLISTER PACKAGE

Cannula is packed in rigid blister and sealed in peel open type tyvek.



Scan to
Visit Product Page



HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.
174, 178/25 BALLABGARH FARIDABAD, HARYANA
121004 INDIA

Manufacturing Site :
HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.
JHARSAINLY-JAJRU ROAD SECTOR-59 PART-II
BALLABGARH FARIDABAD, HARYANA 121004 INDIA

Advena Ltd.
Tower Business Centre, 2nd Flr.,
Tower Street, Swatar, BKR 4013 Malta