


# Cathy<sup>®</sup> NEO

INSTANT-FLSH BK



Getting it Right 1st Time, Every Time

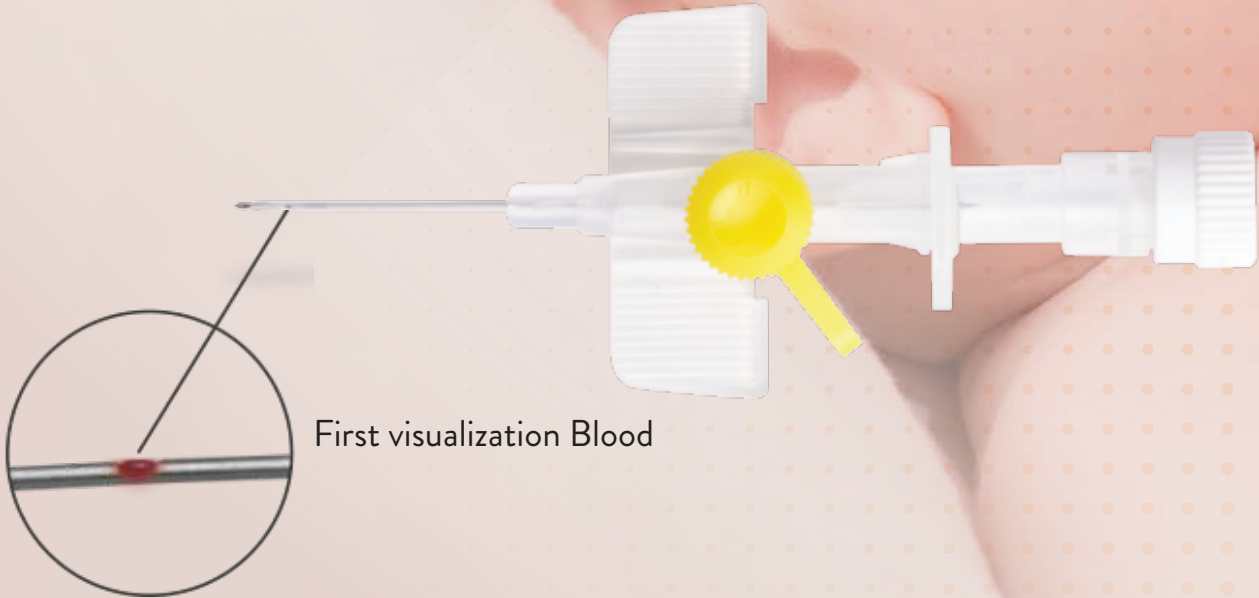

First visualization Blood

WE CARE FOR YOUR LITTLE ONES



# CATHY with “INSTANT FLSH BK”

## Eases Cannulation, Reduces Anxiety

For neonates and children, in 24G & 26G I.V Cannula, a notch has been provided in the needle for “INSTANT FLSH BK” (instant flashback visualisation) at the tip of the catheter. It assists in achieving successful Cannulation in 1st attempt of vein entry. It is specially helpful for neonates with difficult to access veins thus reducing trauma and anxiety – both for you and your patient.

### NEEDLE


Siliconised electrofinished, back cut bevelled stainless steel needle  
Notch on 24G and 26G Cannula

### LUER LOCK CAP

Closure to wing housing permits patient mobility when infusion temporarily discontinued.

### BLOOD STOPPER

The stopper has fine air vent to expedite blood flash back for enabling correct cannulation.

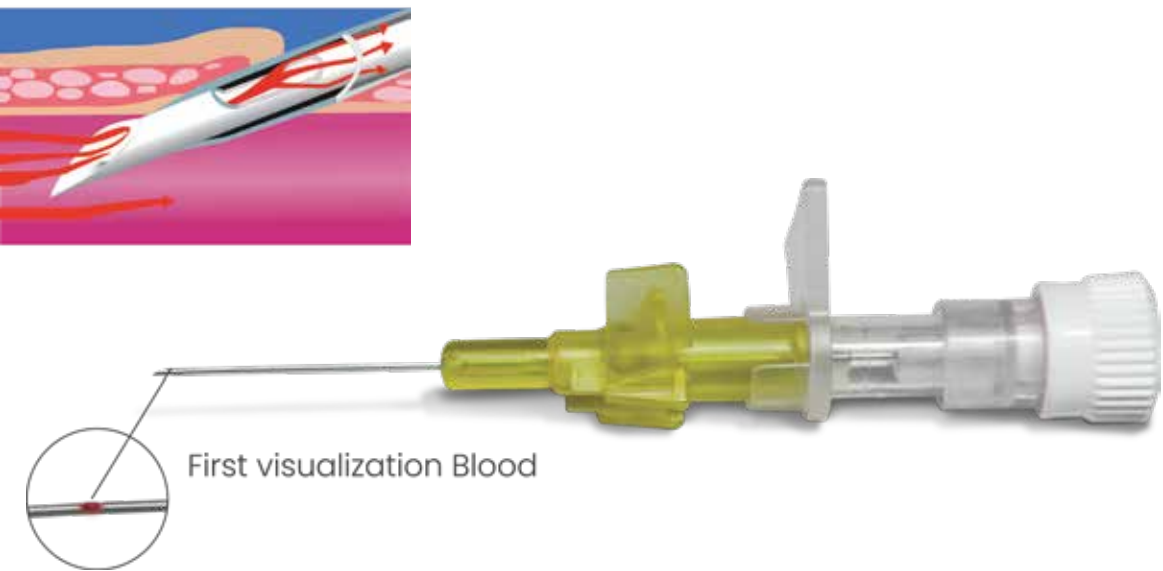

### NEEDLE HUB

Allows blood flash back visualisation.  
Thumb rest for proper grip during insertion.

### CATHETER

Available in VOLEX™ Polyurethane radio detectable 4 strips.

### INJECTION PORT CAP

It facilitates extra medication and also prevents its back flow.

### BLISTER PACKAGE

Cannula is packed in rigid blister and sealed in peel open type medical grade tyvek.

Colour Code	Gauge	Catheter Ext. Dia x Length(mm)	Water flow-rate (ml/min)
Yellow	24G	0.7 x 19	22
Violet	26G	0.6 x 19	15



Scan to Visit Product Video  
Scan to Visit Product Page



**HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.**  
174, 178/25 BALLABGARH FARIDABAD, HARYANA  
121004 INDIA

Manufacturing Site :  
**HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.**  
JHARSAINLY-JAJRU ROAD SECTOR-59 PART-II  
BALLABGARH FARIDABAD, HARYANA 121004 INDIA

www.hmdhealthcare.com | info@hmdhealthcare.com

**Advena Ltd.**  
Tower Business Centre, 2nd Flr.,  
Tower Street, Swatar, BKR 4013 Malta

Revision No:(Jan.2025/Rev.-01)