

VAKU-8™ VAKU-8#

Blood Collection System



With safety shield, offer the ideal solution for accessing difficult veins for blood sampling by utilizing the wings to control needle placement.


The wings ensure no trauma due to the rolling of the needle and minimize chances of vein perforation. Distal vein wall puncture aside from causing patient discomfort due to repeated pricks, can result in increased chances of Phlebitis and other related problems.

A non latex rubber sheath on multiple sample adaptor needle retracts automatically as it penetrates the vacuum tube stopper and recover upon withdrawal. This sheath helps in avoiding accidental needle sticks injury to personnel using it, as well as avoiding blood spillage.

There is a safety shield for covering tube stopper the needle post use to prevent accidental needle stick injuries.

Features

VAKU-8 [‡]	VAKU-8 TM
FLEXIBLE WINGS: Their textured surface easily pinch together to provide a secure grip.	FLEXIBLE WINGS: Their textured surface easily pinch together to provide a secure grip.
NEEDLE: With short bevel, minimises chances of vein perforation.	NEEDLE: With short bevel, minimises chances of vein perforation.
KINK-RESISTANT TUBING: Assures continuous flow.	KINK-RESISTANT TUBING: Assures continuous flow.
FEMALE LUER LOCK: For secure fitment.	FEMALE LUER LOCK: For secure fitment.
PROTECTIVE SHIELD: Ensures there are no secondary needle sticks and no reuse of device.	-----
SELF SEALING NON LATEX PROTECTION SHEATH: Automatically covers needle before & after insertion.	SELF SEALING NON LATEX PROTECTION SHEATH: Automatically covers needle before & after insertion.

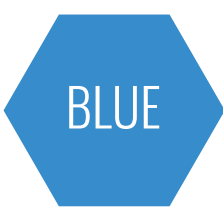

GUAGE	21G	23G	24G	25G
LENGTH	3/4"	3/4"	3/4"	3/4"
PACKAGING	100 Pcs	100 Pcs	100 Pcs	100 Pcs

• VAKU - 8+ IS ALSO AVAILABLE IN 190 MM • OTHER SIZES AVAILABLE ON REQUEST • COLOUR CODE: AS PER ISO: 6009



Scan to
View the product Video
Revision No:(Jun. 2025/Rev-02)



Manufacturing Site:
HINDUSTAN SYRINGES & MEDICAL DEVICES LTD.
99/25 BALLABGARH FARIDABAD, HARYANA
121004 INDIA